



## Data sheet

print date: 2026-01-21

## Fume hood Compact Line DCL-1500 SS WALK-IN



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.


### Standard equipment for fume hoods:

- structure entirely made of galvanized steel without the use of wood-based materials
- toughened glass window, lifted manually with counterweight, opened upwards (max 1850mm). The window has a system preventing the uncontrolled fall of the window
- energy-saving linear LED lighting separated from the fume hood's working chamber with a glass pane
- fume hood main switch
- a double rear wall system with exhaust gaps in the upper and lower part ensuring the correct extraction of vapors that are lighter or heavier than air, without swirls and dead zones
- safety flap in the ceiling ensuring decompression in the event of uncontrolled pressure increase in the working chamber
- full monitoring set (Schneider brand) of air flow velocity in the fume cupboard window compliant with the PE-EN 14175 standard. : PN-EN 14175; PN-EN 14056; PN-EN 13150



## TECHNICAL DATA - DIMENSIONS

overall dims [mm]	kat
width A	1500
height B	2385-2850
depth C	1200
working space dims (mm)	kat
width D	1170
window width D1	1110
height E	2145
depth F	845

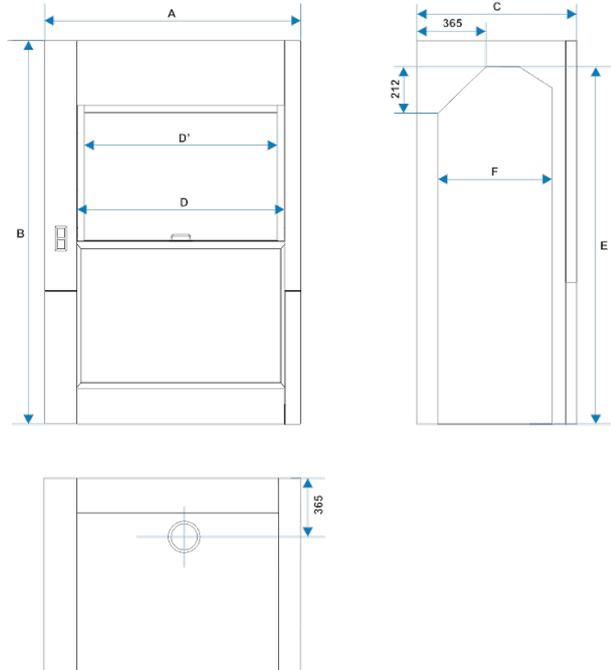

## TECHNICAL DATA

recommended airflow [m3/h]	750 - 1200
required capacity of the ventilation system [m/s]	0,3 - 0,5
nominal power [W]	82
voltage	230 VAC 50Hz
diameter of the air outlet stub [mm]	Ø 250
electrical insulation	class I
standard ventilation system	double rear wall
monitoring set	<i>standard:</i> - FM 550-A-0-E controller with a speed sensor AFS 100 in the fume hood window and the FA-0025-3 operator panel <i>option:</i> - ICM 500-F-EU controller with flow control - ICM 500-FP-EU controller with flow control and cooperation with BMS - other Schneider controllers available on request
raising the window	<i>standard:</i> - manual with counterweight (lock at a height of 500mm, maximum opening 1850mm) <i>option:</i> - SC 500 - Auto Protect system - with foot switch - SC 500 - Auto Protect system - with panel switch
frame and housing material	The frame is made entirely of galvanized steel. External housing made of galvanized and painted with epoxy paint steel
working chamber wall finishing	galvanized steel painted with epoxy paint
warranty	24 months
manufacturer	POL-EKO®

## OPTIONS AND ACCESSORIES



Order number: D-S-EL400

Additional electric socket (400V)





Order number: D-S-EL

Additional electric socket (230V)



Order number: D-WRP

Residual current device



Order number: D-GRZ

Emergency stop button



Order number: D-L-EX-XX.XX

EX Lamp



Order number: D-S-WA

Additional water instalation



Order number: D-S-GA-XXX

Additional LPG instalation



Order number: ICM500-F-EU-XX.XX

Advanced fume hood controller ICM500-F-EU



Order number: ICM500-FP-EU-XX.XX

Advanced fume hood controller ICM500-FP-EU



Order number: D-KRX-XX.XX

Scaffold



Order number: D-SINK-NS

Additional sink



Order number: DCL-LM

Phenolic resin interior



Order number: DCL-PP

Polypropylene interior



Order number: DCL-CR

Buchtal ceramic interior