



Data sheet

print date: 2026-01-21

Fume hood Compact Line DCL-1500 SS



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Compliance with PN-EN 14175 Part 1, PN-EN 14175 Part 2, PN-EN 14175 Part 3, certified by an independent, accredited testing laboratory.

Standard equipment for fume hoods:

- structure entirely made of galvanized steel without the use of wood-based materials
- worktop at a height of 900mm made of monolithic ceramics with a raised rim with a suspended ceramic sink 250x95mm
- toughened glass window, lifted manually with counterweight, opened upwards with a lock at the height of 500 mm (max 810mm), and additionally sliding glass (right / left). The window has a system preventing the uncontrolled drop of the window
- 2 electric sockets 230V 50Hz
- 2 water spouts with valves on the BROEN front panel
- fume hood main switch
- energy-saving linear LED lighting separated from the fume hood with a glass pane
- double system the back wall with exhaust slots in the upper and lower part ensuring proper extraction of vapors lighter or heavier than air without turbulence and dead zones
- safety flap in the ceiling ensuring decompression in the event of uncontrolled pressure increase in the working chamber
- full monitoring set (Schneider brand) air flow velocity in the fume cupboard window compliant with the PE-EN 14175 standard
- CE Declaration of Conformity: PN-EN 14175; PN-EN 14056; PN-EN 13150




POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A.Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

TECHNICAL DATA - DIMENSIONS

overall dims [mm]	kat
width A	1580
height B	2325-2600
depth C	960
working space dims (mm)	kat
width D	1450
window width D1	1265
height E	1220
depth F	635

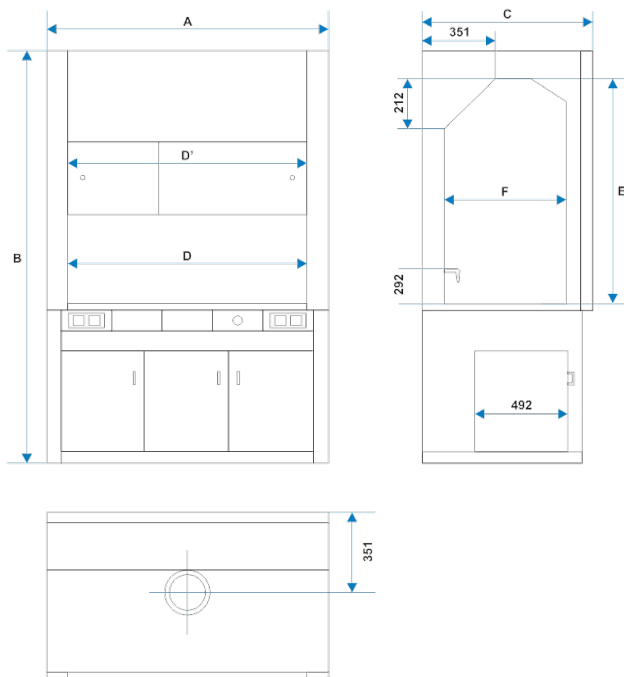

TECHNICAL DATA

recommended airflow [m3/h]	750 - 1200
required capacity of the ventilation system [m/s]	0,3 - 0,5
nominal power [W]	82
voltage	230 VAC 50Hz
diameter of the air outlet stub [mm]	Ø 200
water supply	G 1/2 " male thread
sewage connection diameter [mm]	Ø 50
electrical insulation	class I
standard ventilation system	double rear wall
illumination installed outside the working chamber	LED, class A ++, through insulating glass / from operator panel
monitoring set	<i>standard:</i> - FM 550-A-0-E controller with a speed sensor AFS 100 in the fume hood window and the FA-0025-3 operator panel <i>option:</i> - ICM 500-F-EU controller with flow control - ICM 500-FP-EU controller with flow control and cooperation with BMS - other Schneider controllers available on request
raising the window	<i>standard:</i> - manual with counterweight (lock at a height of 500mm, maximum opening 810mm) <i>option:</i> - SC 500 - Auto Protect system - with foot switch - SC 500 - Auto Protect system - with panel switch
frame and housing material	The frame is made entirely of galvanized steel. External housing made of galvanized and painted with epoxy paint steel
working chamber wall finishing	galvanized steel painted with epoxy paint
worktop	Monolithic ceramics with marine edge
worktop OPTION (available materials)	- Phenolic resin, - epoxy resin - Stainless steel AISI 304, 316
underbench cabinet OPTION (different types)	- metal, ventilated cabinet painted with epoxy paint - cabinet for acids and alkalis made of foamed PVC - explosion-proof cabinets Ascos 30 and 90 minutes
warranty	24 months
manufacturer	POL-EKO®



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

OPTIONS AND ACCESSORIES



Order number: D-S-EL400

Additional electric socket (400V)



Order number: D-S-EL

Additional electric socket (230V)



Order number: D-WRP

Residual current device



Order number: D-GRZ

Emergency stop button



Order number: D-L-EX-XX.XX

EX Lamp



Order number: D-S-WA

Additional water instalation



Order number: D-S-GA-XXX

Additional LPG instalation



Order number: ICM500-F-EU-XX.XX

Advanced fume hood controller ICM500-F-EU



Order number: ICM500-FP-EU-XX.XX

Advanced fume hood controller ICM500-FP-EU



Order number: D-SZ-SS-Z-XX.XX

Underbench cabinet (steel)



Order number: D-SZ-PP-K-Z-XX.XX

Underbench cabinet (polypropylene)



Order number: D-SZ-A90.50-XX.XX

Underbench cabinet (fireproof)



Order number: SC500

Automatic sash window



Order number: D-S-P

Side walls glazing



Order number: D-KRX-XX.XX

Scaffold



Order number: D-SZ-KD-XX.XX

Cuvette for cabinets



Order number: D-SINK-NS

Additional sink



Order number: D-WT-NS

Non-standard worktop



Order number: D-S.C.A.T.

Funnel for gathering leachate



Order number: SC500-L

Automatic sash window with IR sensor