



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

Data sheet

print date: 2026-01-20


Cooled Incubator (ST) ST 1200 Smart PRO



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Advantages of the SMART PRO controller:

- large (7"), clear, full colour touch screen
- LAN, USB ports and WiFi for communication and data transfer
- multi-segment time and temperature programs
- overview of data in tabular and graphic form
- visual and sound alarm
- Admin function for management
- password protected log-in
- internal memory for programs and data storage
- event registry with user notifications
- LabDesk software and user manual for direct download

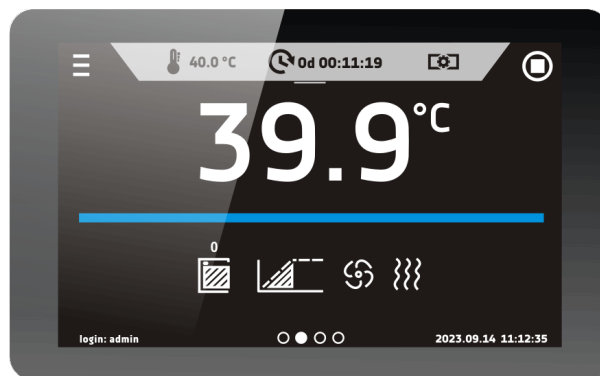

Smart PRO - preview screen



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

TECHNICAL DATA

air convection	forced
chamber capacity [l]	1397
working capacity [l]	1318
controller	microprocessor PID
display	7" full colour touch screen

TEMPERATURE


temperature range [°C]	+3...+70
temperature resolution every ... [°C]	0,1
temperature fluctuation at 37°C [+/°C]*	0,3
temperature variation at 37°C [+/°C]*	1
temperature protection	class 3.3 to DIN 12880

CHAMBER

door type	solid / glass or double (option) /4/
interior	
C Smart PRO	stainless steel to DIN 1.4016
CS Smart PRO	stainless steel to DIN 1.4016
P Smart PRO	acid-proof stainless steel to DIN 1.4301
PS Smart PRO	acid-proof stainless steel to DIN 1.4301
housing	
C Smart PRO	powder coated sheet
CS Smart PRO	stainless steel polished
P Smart PRO	powder coated sheet
PS Smart PRO	stainless steel polished

overall dims [mm] /1/

width A	1470
height B	1990
depth C	900
internal dims [mm]	
width D	1270
width D'	1330
height E	1510
depth F	680
height I	1380

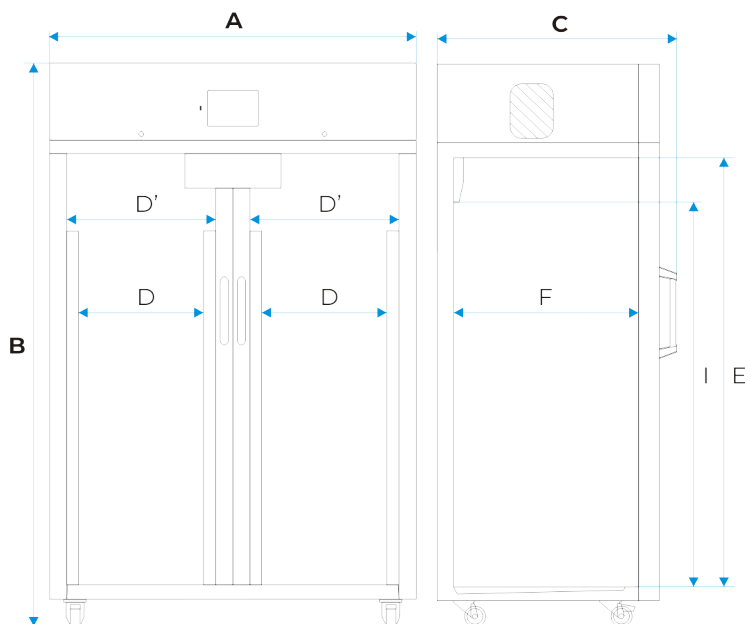

shelves (standard max)	2 x 3 11
max shelf workload [kg] /2/	30
- reinforced shelf version (PW) [kg] /3/	100
max unit workload [kg]	300
weight [kg]	221



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

ELECTRICAL PARAMETERS

voltage	230V 50/60Hz
nominal power [W]	650
refrigerant	R290 / GWP=3
warranty	24 months
manufacturer	POL-EKO®

all the above technical data refer to standard units (without optional accessories)

* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

$K = \pm (T \text{ average max.} - T \text{ average min.}) / 2$

** - other power supplies on request


1 - depth doesn't include 50 mm of power cable, the width does not include the 20 mm of rubber plug

2 - on uniformly loaded surface

3 - reinforced shelf


4 - additional internal glass door

OPTIONS AND ACCESSORIES

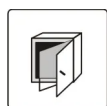

Order number: */C

Internal glass door

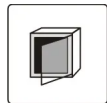


Order number: */A

External glass door




Order number: */P INOX

Stainless steel wire shelf INOX

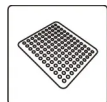


Order number: */PP

Perforated shelf

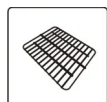


Order number: */PW

Reinforced shelf

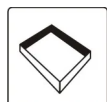


Order number: KUV GN*/*

Stainless steel cuvettes

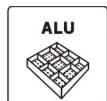

Order number: ST/CHL/SWP ALU

Aluminum drawer with powder coated slides

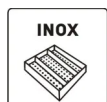


Order number: ST/CHL/SWP INOX

Stainless steel drawer with powder coated slides



Order number: ST/CHL/SWP INOX

Stainless steel drawer with stainless steel slides

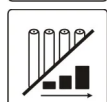

Order number: */FIT


Phytotron system



POL-EKO®


(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

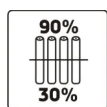


Order number: FIT/R3

FIT panels independent control

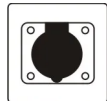


Order number: GNZ

Internal socket

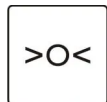

Order number: BRT/*/L or IQ/OQ/PQ

Calibration and IQ, OQ, PQ qualification

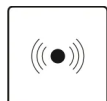

Order number: RFID LOCK (SMART PRO)

RFID LOCK (SMART PRO)




Order number: BPP 12

Battery backup for display




Order number: PORT ALARM

Dry alarm contact

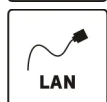

Order number: LANK

LAN cable

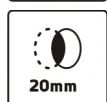

Order number: OCZ/20

Non-standard access port 20 mm



Order number: OCZ/30

Additional access port 30 mm



Order number: OCZ/60

Non-standard access port 60 mm



Order number: OCZ/100

Non-standard access port 100 mm